## 18" Mini Portable Clean Room Hood

## Model # SS-118-IPURE

The 18" ISO-Pure Mini Portable Clean Room Hood provides a solution for applications that require a low cost, compact, and clean positive-pressure environment. This system offers up to an ISO Class 4 clean room cabinet to protect the process from contamination. It utilizes a dual stage filtration process, consisting of a pre-filter and main filter (choice of **HEPA** or **ULPA filter**), to remove particulate from the ambient air and provide filtered air to the workspace. The **HEPA filter** provides up to 99.97% efficiency on particles down to 0.3 microns and larger, while the **ULPA filter** offers up to 99.9995% efficiency on particles down to 0.12 microns and larger. This mini portable clean room hood, while not intended for sterile compounding, has been utilized in variety of applications including stem cell therapy practices, medical applications, tissue culturing, data recovery,

This unique system features a small footprint, energy efficiency, and low maintenance operation. The **18**" **ISO-Pure Portable Clean Room Hood** can be easily integrated into existing workspaces because of its size and can be transported to various locations without difficulty because of its lightweight and portable nature.

cleaning processes, manufacturing, and much more.



MADE IN

THE USA

SINCE 1983

18" ISO-Pure Mini Portable Clean Room



18" ISO-Pure Mini Portable Clean Room can be a compact solution for data restoration.

## **Product Features**

- Provides up to an ISO Class 4 Clean Room
- · Includes: anti-static vinyl curtains
- Highly portable and lightweight
- Quiet operation
- Reliable, low maintenance operation
- · Long filter life
- Simple "no tool" filter change



21221 FM 529 Rd, Cypress, TX 77433 1.800.799.4609 • sales@sentryair.com

**Outer:** 19" x 20" x 24"

Interior: 18" x 18" x 15.5"

**Access Area:** 7.5" x 16"


**Hood Dimensions** 

 $(D \times W \times H)$ 

Warranty



Dimensions are approximate

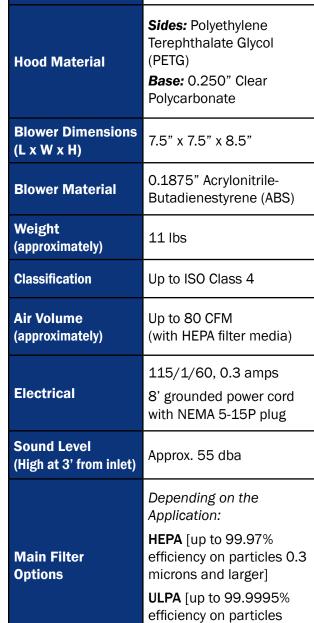

0.12 microns and larger]

Limited three-year

to materials or

workmanship.

warranty from date of

shipment on defects due



21221 FM 529 Rd Cypress, TX 77433 1.800.799.4609 sales@sentryair.com

www.sentryair.com