

24" Portable Clean Room Hood

Model # SS-324-PCR

The **24" Wide Portable Clean Room Hood** offers up to an *ISO Class 5 cleanroom* to help protect the process from outside contamination and particulate. This system utilizes a vertical laminar airflow pattern and a dual-stage filtration to create positive-pressure environment by supplying filtered air to the workspace. The **Portable Clean Room Hood** draws in ambient air and processes the air through a pre-filter and main filter. Main filter choices include **HEPA filtration** (up to 99.97% efficiency on particles down to 0.3 microns) and **ULPA filtration** (up to 99.9995% efficiency on particles down to 0.12 microns).

The **Portable Clean Room Hood** offers a reliable and low maintenance operation for a variety of applications including *IV infusion, parts inspection, optics/laser assembly, tissue culturing, data recovery, stem cell therapy practices and other medical applications*. This unit's compact size and portable nature also make it a good choice for field-testing applications that require frequent travel. To manage the air volume, this unit comes standard with a variable speed controller. Also included is a LED light with an optional UVC light for surface disinfection.

Our line of portable clean rooms is not intended for chemotherapy.

Product Features


- **Includes:** anti-static vinyl curtains, LED light, and variable speed controller (115V only)
- **Optional:** UVC light, Magnehelic Gauge, hour counter, and removable spill tray
- **Standard Sizes:** 12", 18", 24", 30", and 40"
- Provides up to an ISO Class 5 Clean Room
- Highly portable and lightweight
- Quiet operation
- Reliable, low maintenance operation
- Long filter life
- Simple "no tool" filter change



24" Portable Clean Room Hood



24" PCR provides a clean workspace for tissue culturing.

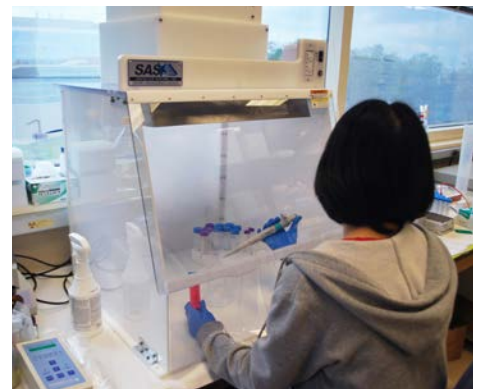

The 24" PCR Hood provides a sterile workspace for a variety of medical applications



21221 FM 529 Rd, Cypress, TX 77433
1.800.799.4609 • sales@sentryair.com

www.sentryair.com

Positive-Pressure Clean Room

Product Specifications



Dimensions are approximate

Hood Dimensions (D x W x H)	Outer: 24" x 26" x 25" Interior: 23.5" x 24" x 24" Access Area: 8" x 24"
Hood Material	Sides: 0.25" white acrylic Base: 0.5" white high density polyethylene Top: 0.5" clear polycarbonate Front Window: 0.25" clear polycarbonate
Blower Dimensions (L x W x H)	15.5" x 13" x 14.5"
Blower Material	16 ga. carbon steel with polyurethane epoxy coating on all surfaces
Weight (approximately)	95 lbs
ISO Classification	Up to ISO Class 5
Air Volume (approximately)	Up to 350 CFM (with HEPA filter media)
Electrical	115/1/60, 2 amps 220/1/50, 0.75 amps 8' grounded power cord with NEMA 5-15P plug
Sound Level (High at 3' from inlet)	Approx. 55 dba
Pre-Filter	MERV 8 Pre-Filter
Main Filter Options	<i>Depending on the Application:</i> HEPA [up to 99.97% efficiency on particles 0.3 microns and larger] ULPA [up to 99.9995% efficiency on particles 0.12 microns and larger]
Warranty	Limited three-year warranty from date of shipment on defects due to materials or workmanship.
Patent	#5,843,197



21221 FM 529 Rd
Cypress, TX 77433
1.800.799.4609
sales@sentryair.com

www.sentryair.com