

# Laboratory Equipment

## Fume Hood Base Cabinets

Base cabinets are available to fit under standard laboratory work surfaces and counters, with safety features suitable for hazards as specified.


**Standard** cabinets follow the same design and construction as 7000 series cabinets.

**Acid** cabinets follow the same design and construction as 7000 series cabinets but include an interior lining of Wesliner material and front louvers for venting. "ACID" label is included.

(See page 39 for more info)


**Flammable Storage** cabinets feature double wall construction and rear vents. Available with Manual or Self-close (Sliding) doors. OSHA and NFPA compliant. (Units are finished in SAFETY YELLOW unless otherwise specified)

**Standard**

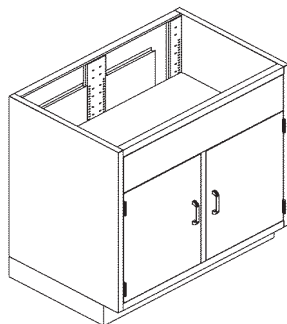


Model	Dim. W x H
7119-29FH	29" x 35-3/4"
7119-35FH	35" x 35-3/4"
7119-47FH	47" x 35-3/4"
7119-58FH	58" x 35-3/4"

**Acid Storage**

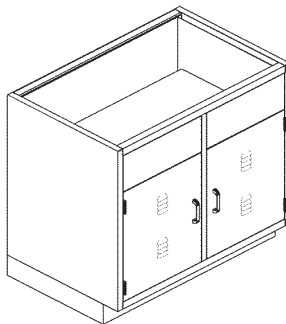

Model	Dim. W x H
7000-29FH	29" x 35-3/4"
7000-35FH	35" x 35-3/4"
7000-47FH	47" x 35-3/4"
7000-58FH	58" x 35-3/4"

**Flammable Storage**

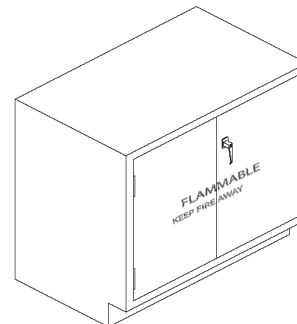

Model	Dim. W x H x D
L124 Manual	35" x 35-9/16" x 22"
L224 Self-Close	35" x 35-5/8" x 22"
L136 Manual	47" x 35-9/16" x 22"
L236 Self-Close	47" x 35-5/8" x 22"
L144 Manual	59" x 35-9/16" x 22"
L244 Self-Close	59" x 35-5/8" x 22"

### Vent Kits

#### **Model VK-A**

Vent Kit for 7000 series Acid cabinets. PVC Vent extends from the cabinet to the worktop so fumes can be exhausted by the hood. Kit includes elbow and two threaded connections.

#### **Model VK-F**

Vent Kit for Flammable storage cabinet. Flexible plastic hose extends from the cabinet directly into ductwork, bypassing the fume hood entirely so dangerous fumes are exhausted directly into the ductwork. Kit includes 10' hose and adapters.

\*Other materials available upon request.

